

WHEN YOU PUT PATIENTS' HEALTH FIRST,  
REMARKABLE THINGS FOLLOW



State-of-the-art manufacturing facilities.



### OUR RANGE OF PRODUCTS INCLUDE

SINGLE USE CIRCULAR  
STAPLER FOR RECTAL  
PROLAPSE AND  
HEMORRHOIDS

SINGLE USE  
RELOADABLE LINEAR  
STAPLER

SINGLE USE CURVED  
INTRALUMINAL  
CIRCULAR STAPLER



CHEX<sup>™</sup>  
HEALTH FIRST

INFORMED DESIGN.  
SCRUPULOUSLY APPLIED.  
METICULOUSLY PRODUCED.

## SINGLE USE RELOADABLE LINEAR CUTTER STAPLER



INFORMED DESIGN.  
SCRUPULOUSLY APPLIED.  
METICULOUSLY PRODUCED.



FRANKENMAN INTERNATIONAL LTD  
SUITE B 13/F WING TAT COMMERCIAL BUILDING  
121-125 WING LOK STREET SHEUNG WAN HONG KONG  
T +852 3106 3035 F +852 3585 0519  
E info@frankenman.hk

[www.frankenman.hk](http://www.frankenman.hk)

Chex and Chex Healthcare are registered trademarks  
of Frankenman International Limited

CHEX<sup>™</sup>  
HEALTH FIRST



A RANGE OF PRODUCTS TO MEET THE NEEDS OF MODERN SURGICAL PRACTICE

The Chex™ range of Linear Cutter Staplers focuses on reliability and functionality. Designed to provide accurate and reliable staple formation and tissue resection... time and time again.

At Frankenman, we believe we're not so much in the business of making instruments, but in the business of making surgeries more successful.

That's why we listen to doctors who've put our instruments to use. We read their reports. We invite them to think with us, to brainstorm, to help us develop new and innovative ideas.

That way we know we will be able to deliver the highest quality products to the medical profession both today and way into the future.



THE CHEX™ LC RANGE  
ACTUAL SIZE



SIDE VIEW



KEY FEATURES

- 1  
CLEAR GRADUATIONS AND INDICATION OF 'CUTTING ZONE' ON BOTH SIDES OF STAPLE HOUSING.
- 2  
TISSUE RETAINING PIN.
- 3  
UNIQUE 90 DEGREE RELOAD 'SHOULDER' TO ACT AS A POSITIVE TISSUE LIMITER TO AVOID CONCERTINA EFFECT OVER RELOAD AND RISK OF PRE-CUT INTERNALLY.
- 4  
ERGONOMIC AND RIDGED HANDLES TO PREVENT SLIPPING DURING FIRING.
- 5  
NEWLY DESIGNED FIRING SLIDER.

- 6  
SPECIALLY DESIGNED WIDER AND MORE OPEN JAW MECHANISM.
- 7  
NEW SAFETY GATE MECHANISM TO PREVENT FIRING WITH EMPTY CARTRIDGE (PAT. PENDING).

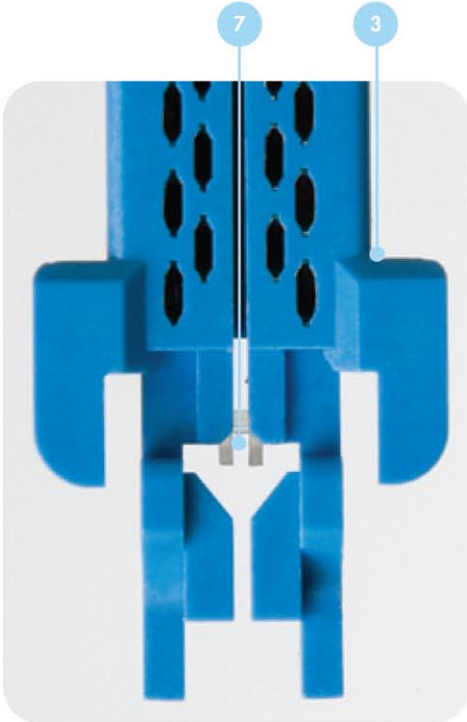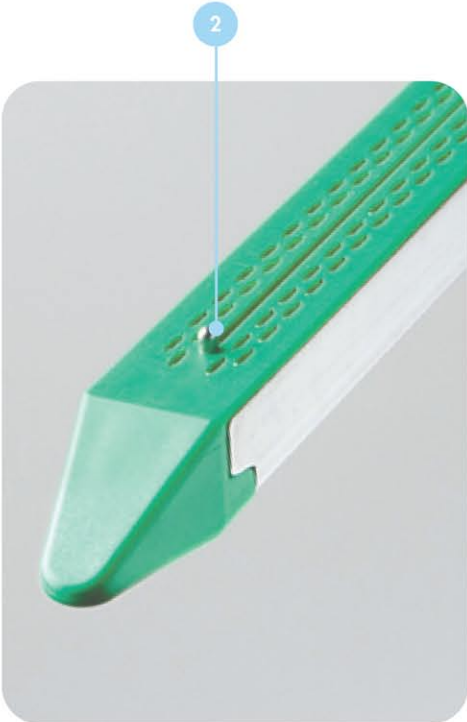

PRODUCT SPECIFICATIONS

PRODUCT CODE				CLOSED STAPLE SIZE	STAPLING LINE LENGTH	STAPLE QUANTITY
LINEAR CUTTER STAPLER	CARTRIDGE RELOADS					
LC60G	GREEN	LCR60G	GREEN	2.0MM	61MM	60
LC60B	BLUE	LCR60B	BLUE	1.5MM	61MM	60
LC80G	GREEN	LCR80G	GREEN	2.0MM	81MM	80
LC80B	BLUE	LCR80B	BLUE	1.5MM	81MM	80
LC100G	GREEN	LCR100G	GREEN	2.0MM	100MM	100
LC100B	BLUE	LCR100B	BLUE	1.5MM	100MM	100

KEY FEATURES

- CAM GEAR DESIGN: ENSURES CONSISTENT STAPLE FORMATION HELPING PREVENT BLEEDING.
- EMPTY CARTRIDGE SAFETY LOCK: CUTTING BLADE WILL ONLY ACTIVATE WITH FULL CARTRIDGE.
- RELOADABLE COLOUR CODED CARTRIDGES.